BEAUTIFLY

Luminare Rejuva Skin Visage

Golden

series

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- www.beautify.eu
- www.youtube.com/beautifyeu

Rated voltage:	AC100-240V
Rated frequency:	50/60Hz
Rated power consumption:	120W
Dimensions:	Dimensions ~ 325mm x 82mm x 101mm (when the handle is folded)
Libra:	Weight ~ 1240g
Power cord:	1.5m
Handle connection cable:	1.3 m
Timer:	Automatic shutdown after 20 minutes
Use environment:	Temperature 5°C - 35°C, relative humidity ≤75%
Material:	Main case: PC + ABS, silicone, aluminum alloy; Various attachments: PC + ABS, aluminum, glass

Description of parts - device in the base [see page 2, f2]

- Level LEDs - The skin color sensor detects 10 levels and automatically adjusts the flash power level. All LEDs are always on when the device is placed on the base.
- Status LED - You can check what mode is set by the color of the LED. The LED is always on in standby mode. When it is placed on the skin, it will flash.

Modes	Color of the lamp	Blink level
Sensitive skin	Blue	Short
Normal skin	White	Mediocre
High power	Orange	Tall

- Power/Mode button; Long press: Power on/off. Short press: Each time the button is pressed, the mode is switched. Sensitive skin mode, normal mode and high power mode can be selected.
- Basis
- Connecting cable.

Description of parts - device disconnected from the base [see page 2, f3]

- Air outlet
- Touch sensor
- Cooling surface
- Flash outlet
- Skin color sensor
- Flash button - Short press: single flash, Long press and hold the button: flashes for about 3 seconds.
- Air inlet

Description of parts - stand and cable [see page 3, f5]

- UV light - When the power cord is plugged into the socket or the device is suspended, the stand enters standby mode and the UV light works automatically for about 3 minutes. When the UV light is working, all level lights flash slowly.
- Stand
- Power cable
- EU power plug

Accessories [see page 3, f6]

- Face/finger cap
- VIO overlay
- DPL skin care pad
- NIR overlay
- Safety glasses

How to install/remove the overlay [see page 3, f7]

The cap is of the magnetic type.

Install the cap so that it covers the flash exit surface. When disconnecting, hold both sides of the cap and remove it gently

IPL function (Intense Pulsed Light - VIO overlay can be used to care for the VIO zone; face/finger overlay)

Hair removal is performed using the IPL function. The pulse intensity level can be selected from three modes (sensitive skin, normal skin, high power mode) according to the skin condition, color and area to be used.

DPL function (Dye Pulsed Light - DPL care pad)

Use for broken capillaries (redness) and enlarged pores (shrink pores); Pigmentation Reduction: DPL can target and reduce the appearance of pigmentation problems such as age spots, sun spots and freckles. It works by breaking down excess melanin in the skin, leading to a more even and balanced skin tone. Skin Firming: DPL technology can help to firm and tighten the skin by stimulating the production of collagen and elastin, which helps improve skin elasticity and reduce sagging. Skin Renewal: DPL can help stimulate the production of collagen, which is essential for maintaining youthful and healthy-looking skin. Targets fine lines, wrinkles and uneven skin tone, leaving you with a more radiant and refreshed complexion. Relief from redness and rosacea: DPL light energy can help calm and reduce redness associated with conditions such as rosacea. It works on the blood vessels responsible for redness, leading to a more even skin tone and reduced redness over time. How often to use DPL: no more than twice a week.

NIR function (Near InfraRed - NIR care pad)

Skin Firming and Tightening: NIR technology stimulates the production of collagen and elastin, leading to improved skin elasticity and firmness. This helps tighten loose skin, reduce sagging and promote a younger appearance. Wrinkle Reduction: By enhancing collagen synthesis, NIR can help reduce the appearance of fine lines and wrinkles. The heat generated by NIR penetrates deep into the skin, stimulating the production of new collagen fibers, resulting in smoother and more youthful skin.

Skin Rejuvenation: NIR therapy can improve the overall texture and tone of the skin. It targets areas of uneven pigmentation, sun damage and age spots, promoting a more even and radiant complexion. How often to use NIR: no more than once a day, but not together with DPL.

Continuous flash mode

It is recommended to use this mode for large surfaces such as legs and arms.

Hold the flash button for about 3 seconds to switch to continuous flash mode.

Use it by moving it across your skin without removing the device from your skin.

The continuous flash mode will continue to operate as long as it is in contact with the skin.

If the device comes away from your skin, press the flash button again. If skin color and touch sensor are not detected for approximately 3 seconds, continuous flashing will stop.

Safety features

Skin color sensor (skin color recognition)

The skin color sensor reads the skin color of the flash surface in 10 levels and automatically adjusts the flash level. Skin color recognition does not work when using overlays.

Touch sensor

When the flash output is fully in contact with the skin, it is possible to fire. While the status light is flashing.

Cooling function

The cooling surface is located around the flash outlet. Cools the skin during the treatment.

Auto power off function

When the power is turned on and there is no operation for 3 minutes, the power will be turned off automatically. Moreover, after 20 minutes of use, the power will also be turned off automatically. The fan will run for approximately 30 seconds after the power is turned off.

Application points [see page 4, f8]

Place the flash output close to the skin and glide gently over the skin.

- Move away from the skin
- NOTE: Make sure there are no gaps.

Care areas (points)

Hair removal (IPL) [see page 4, f9]

- Armpit / b. Chest / c. Arm / d. Arms and fingers / e. Leg / f. Face / g. VIO zone

Skin care (GLP) [see page 4, f10]

- Face / b. Neck and cleavage

Methods of use [see page 4, f10]

Hair removal (IPL), skin care (DPL)

- Connect the power cable to the base and a wall outlet.
- Press and hold the Power/Mode button [see page 2, f1-a] for about 3 seconds to turn on the device. The status LED will light up and the fan will start operating (in standby mode), [see page 2, f2]:
 - Level LEDs
 - Status LED - on; in standby mode / flashing: flash available
 - Power/Mode button

- Briefly press the Power/Mode button [see page 2, f1-a] to switch mode. You can check the mode by the color of the status light. The flash level increases in the order: sensitive skin mode → normal mode → high power mode.

- Place the flash output on your skin (without the cap), then briefly press the flash button [see page 2, f1-b] when the status LED begins to flash. The skin color sensor automatically identifies skin color in 10 levels and adjusts the appropriate level/diode. To activate continuous mode, press and hold the flash button [see page 2, f1-b] for about 3 seconds.

Single flash (shot) - press the flash button once for a given position. Then change position and "shoot again"

Continuous flash - when caring for wider surfaces, such as arms and legs, move the device over the skin.

The continuous flash mode lasts as long as it is in contact with the skin. If you briefly press the flashing button [see page 2, f1-b] again or leave it without skin contact for approximately 3 seconds or longer, the continuous flash mode will turn off.

- After the procedure, turn off the device - press and hold the power/mode button [see page 2, f1-a] for about 3 seconds. Hang the device on the base (the device will enter standby mode). The UV lamp automatically turns on for approximately 3 minutes. After 20 minutes of use, the power will turn off automatically.

- Make sure the device's fan has stopped and unplug the power cable from the base and power socket. After care:
 - After care, moisturize the skin thoroughly.
 - If the skin is hot after treatment, cool it down with an ice pack or a cold towel.
 - After care, avoid exposing your skin to direct sunlight and be careful not to get sunburned during care.
 - After each use, remove sebum stains or the cosmetic fluid included with the product for maintenance purposes.
 - Before maintenance, remove the plug from the socket.
 - After each use, wipe the flash outlet and skin color sensor with a dry cloth.

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For the neck area, place the treatment window in a horizontal position and move the device from top to bottom. For the décolletage area, place the treatment window in a vertical position and move the device from the inside to the outside.

1. Take care of both sides of the skin sequentially from 1 to 6, moving continuously, leaving no gaps.
2. Use on bare skin (after removing makeup) on your face.
3. Continuous flash mode can be used on the face, but it is recommended to use single exposure mode when caring for the skin around the eyes. Close your eyes when exposing your face.
4. Do not use on eyebrows, eyelids or lips. When an attachment is attached, use it at an exposure level where you can feel a gentle warmth.
5. When using the device with the overlay, you may feel a slight warmth.
6. After the treatment, moisturize the skin thoroughly.

Warnings

- Please remove any metal accessories (earrings, rings... ect) before use.
- The flash is very bright and may damage your eyes
- Do not direct the flash into your eyes
- Do not look directly at the flash part
- Do not use if damaged or broken.
- Do not use repeatedly on the same area more than twice in one day
- Prolonged or repeated flashes in the same place may cause skin damage
- Do not use with blocked air inlet/outlet
- Leave the device for 15 minutes or more before using it again
- Do not use the device on clothes, hair or other objects
- Use the device at temperatures between 5°C and 35°C
- If the device has not been used for a long time, check that it is working properly before using it again
- Do not use power cables other than the dedicated ones
- Do not leave the device switched on
- Do not allow small items, foreign objects, dirt, moisture or dust to stick to the main unit, gaps, ventilation holes, etc.
- Do not disassemble, repair or modify the device
- If the device or accessories emit an unpleasant odor or heat, please stop using them
- If any abnormality occurs, immediately turn off the power of the main unit and remove the plug from the socket
- Do not hit the device against other objects or things
- If the device is dropped or damaged, turn off the power and unplug the device
- Do not pull or coil the power cord
- Always turn off the power and unplug the power plug when you are not using the device.
- When removing the plug from the socket, hold both sides of the plug.
- When using on the face, make sure to remove any makeup, sunscreen, perfume, etc. and use on exposed skin.
- Avoid use immediately after sunbathing. Do not use together with a peeling or care mask.
- When using for the first time, please familiarize yourself with the device, use it for a short time at low level (sensitive skin mode). Even if you do not feel pain, it does not increase strength and use the device at the appropriate level, then gradually adjust the level
- For 24 hours after treatment, be careful not to expose the treatment area to ultraviolet (UV) light. If redness is still present after 24 hours, avoid UV radiation. Also use sunscreen during the treatment.
- After each use, clean the device and remove hair and dirt.

It is prohibited to use the device by:

- Kids
- Menstruating, pregnant and breastfeeding women
- People suffering from heart disease or other cardiovascular conditions
- People suffering from epilepsy, autoimmune diseases, and photosensitivity
- People with cancer, hematopoietic disorders, circulatory disorders, severe inflammatory symptoms, abnormal blood pressure, blood clots
- People with skin abnormalities such as wounds, eczema and swelling, atopic dermatitis and people with very sensitive skin.
- People who use dermatological treatment for acne or other skin ailments
- People with subfebrile or feverish condition
- People who do not feel temperature or external stimuli.
- People with physical disabilities
- People taking medications
- People under the influence of alcohol
- People with wet or damp hands

It is prohibited to use the device in areas such as:

- Lips, pubic area and mucous membrane, tanned skin, eyebrows, hair, eyeball, eyelids, eye area, throat, ear, nipple, areola, navel area, tattoo
- Areas where blood vessels are visible, such as the vena cava
- Scars, melasma, discoloration and birthmarks
- Places after aesthetic medicine treatment (about 2 weeks after surgery)
- Wet or moist skin

Storage:

- Keep the device away from children
- Keep the device away from water, humid places, steam, dusty places, exposed to direct sunlight
- Keep the device away from flammable materials (alcohol, thinner, gasoline, nail polish remover, spray, etc.)
- Keep the device away from places where strong radio waves or magnetic fields are generated
- Keep the device away from heating devices

Cleaning

- For daily maintenance, use neutral detergent and a soft cloth. Do not clean the device with organic solvents or other chemicals. Do not rinse with water.
- Avoid high-pressure washing and cleaning with hard objects that may scratch its surface.
- Between treatments on different parts of the face and body, disinfect the tips using agents intended for beauty devices

More information about the company and products [page 2 - f4]:

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
[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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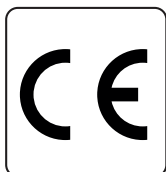[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautifly Sp. z o.o. Plac Bankowy 2, 00-095 Warszawa, Polska
Model: Luminare Rejuva Skin Visage Golden